
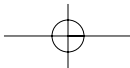
Tizian Blank High translucent

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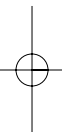If there is anything in this instruction for use that you do not understand, please contact our customer service department before using the product.

As the manufacturer of this medical device, we inform our users and patients that all serious events occurring in connection with it must be reported to us (the manufacturers) as well as the relevant authorities in the Member State where the user and/or patient is resident.

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Instructions

Description / indications

Tizian Blank High translucent is made of a highly translucent zirconium dioxide material and is indicated for the production of full zirconium individual crowns, partial crowns, max. 3-unit bridges (incisors and premolars), inlays, onlays, veneers as well as implant-supported individual crowns on supraconstructions for anterior and posterior teeth. Tizian Blank High Translucent can be designed in a reduced shape to be veneered with ceramics or it can be designed in a fully anatomical fashion.

Safety note

Tizian Blank High Translucent is supplied in a presintered condition. The blanks show only a limited solidity and exhibit a residual porosity. Therefore, it is necessary to handle them carefully. Please check immediately after receiving the consignment whether packaging and product are intact.

Attention: The pre-sintered structure must never come into contact with water or steam.

General information

The information provided in this instruction manual must be passed on to any person using the products mentioned therein. The products must only be used by qualified personnel. The user is obliged to use the products in accordance with the present instruction manual and with appropriate hygiene measures and to verify on his / her own responsibility whether the products are suitable for the individual patient situation. The user will be held fully responsible for the appropriate and correct use of the products. The manufacturer assumes no liability for incorrect results in form of direct or indirect damages or any other

damages that occur from the use and /or the processing of the products. Any claim for damages (including punitive damages), is limited to the commercial value of the products. Independently of this, the user is obliged to report all serious incidents that occur in connection with the products to the competent authority and to the manufacturer.

Contraindications

Should intolerance to the products be suspected, use them only after an allergy screening test has been performed with issue of medical statement of non-allergy conditions.

Processing of Tizian Blank High Translucent

1. Preparation guidelines

Chamfer or shoulder preparations are appropriate for zirconium dioxide. Unsuitable are e. g. channel preparations or tangential preparations. Those preparation types could lead to failure. The following preparation angles are recommended: horizontal 4 - 6° and vertical 3 - 5°. All connecting areas must be designed in a rounded and smooth fashion. Avoid sharp angles. Provide an optimal insertion vector between the posts. Remove at least 1.5 mm of substance from the occlusal surfaces of posterior teeth and 2 mm from the incisal surfaces of anterior teeth.

2. Dimensioning of the framework

Before sintering, the wall thickness must be at least 0.7 mm. Finally, after dense-sintering, the crown walls must not be thinner than 0.4 mm. The margins should be at least 0.2 mm thick. The connector dimensions in the anterior area should not fall short of 3 x 4 mm in the anterior sector and 4 x 4 mm in the posterior sector. The connector geometry is of greatest

significance for the crack resistance of zirconium dioxide frameworks. Therefore, choose the greatest connector height possible.

In order to avoid the so-called „chipping“, design the framework in an anatomical fashion allowing a layer thickness of 1-2 mm for the ceramic veneer.

3. Milling

Only use machines, systems and tools to process Tizian Blank High Translucent which are suited to milling presintered zirconium dioxide materials. The milling machines must have an appropriate mounting fixture for the blanks. The blanks are dry-milled underneath a local extractor system. Please adhere to the instructions of the milling machine's manufacturer.

4. Removing the restorations from the blank/finishing

- Carefully cut the milled restorations from the blank with help of diamond mills.
- Finish the restoration with clean, cross-cut hardmetal mills.
- Avoid thermal damage by ensuring to use low revolutions and minimal pressure, and by cooling sufficiently.
- Blow off the restoration with oil-free compressed air.

5. Coloring

All conventional zirconium dioxide coloring liquids (e.g. Color Liquid by Zirkozahn) can be used. Please refer to the manufacturer's instructions. Pay special attention to the drying times.

Obtain best possible coloring results by individually painting the restoration with a brush. Dipping the restoration into the stain is not recommended. The painting technique is similar to the technique used before glaze firing. Use a brush made only of plastic or wood, without any metal parts, as metal parts could cause discoloration. Please observe the instructions of the stain manufacturer.

6. Sintering in a high temperature furnace

Remove all milling dust from the frameworks before dense-sintering. By sintering, the restoration which was produced about 25% bigger than necessary, will be shrunk to its final size. During sintering, the restoration will obtain its final characteristics. It is absolutely necessary to use a suitable furnace which is approved for presintered zirconium dioxide objects.

Sintering program:

	Weight (tooth unit + sintering bar)		
	0,0 - 2,0 g	2,0 - 4,0 g	> 4,0 g
start	room temperature	room temperature	room temperature
heating rate	max. 14 °C/min.	max. 8 °C/min.	max. 5 °C/min.
final temperature	1450 °C	1450 °C	1450 °C
holding time at final temperature	90 min.	120 min.	120 min.
cooling rate	8 °C/min.	8 °C/min.	5 °C/min.
Open the furnace at	room temperature*	room temperature*	room temperature*

* Higher temperatures than room temperature are possible (max. 300 °C). In isolated cases, this may lead to cracking.

Ensure that wide-span constructions are well supported, thereby avoiding sintering distortions. We recommend removing the restorations from the furnace at room temperature.

7. Finishing the restorations

After sintering, place the restoration on the die, check the margins and, if necessary, finish with a rotating diamond tool under water cooling.

If at all possible, do not finish dense-sintered objects mechanically. If mechanical finishing cannot be avoided, please only use immaculate water-cooled diamond tools to avoid damaging the structure of the object. Avoid extreme local heat (e. g. when using steam jets or when polishing to high lustre).

Do not apply excessive pressure during finishing and take special care when working with very graceful constructions. Ensure that the wall thickness does not fall short of 0.4 mm. Preferably, do not grind the interdental connectors. In principle, do not grind the connectors from basal direction or apply chamfers. Avoid sharp edges and the use of separating discs.

8. Veneer

Only use veneering ceramics suited to veneering zirconium dioxide, e. g. Tizian ceramic (Schütz Dental). Please adhere to the manufacturer's instructions.

To obtain an optimal bond, we suggest a wash firing with dentin mass (100 °C higher than the normal firing temperature).

Zirconium dioxide is a poor heat conductor. In consequence, it takes significantly longer to fire the ceramics on heavy bridge units than on crowns. This fact leads to a special firing procedure:

Weight per unit	Heating rate
0 - 1 g	55 °C/min.
1 - 2 g	45 °C/min.
2 - 3 g	35 °C/min.
3 g	25 °C/min.

Attention: The more solid a construction, the lower a heating rate must be selected. Cool slowly to avoid cracks. Ensure that the minimum holding time of 2 - 3 min. is kept in any case. Remove frameworks from the furnace only after the temperature has fallen below 200 °C. Do not place warm constructions onto a cold surface (e. g. table top), as this can cause cracking, as well.

9. Fixation

Fix the restoration with a translucent fixing composite. If an opaque fixing composite is used, the primary structure may gleam through. Only use materials to integrate the restoration which are compatible with zirconium dioxide. Observe the manufacturer's instructions.

Please note!

When milling partially sintered or dense-sintered zirconium dioxide, respirable dust can be discharged which may damage the lungs or lead to irritations of the eyes and skin. Therefore, zirconium dioxide must only be processed while using a proper and efficient extractor system. Technicians are recommended to wear protective masks and goggles.

Storage

Store the blanks in a dry place at temperatures between 10 °C and 50 °C. Do not subject the blanks to vibrations or thrusts. Generally avoid any contamination.

Disposal

The material can be disposed of with normal household waste.

Guarantee

Our processing instructions, regardless of whether they are issued verbally, in writing or in the form of a practical demonstration, are based upon our own experiences and can therefore be regarded solely as guidelines.

Side-effects

With proper use of this medical device, unwanted side-effects are extremely rare. Reactions of the immune system (e. g. allergies) or local discomfort, however, cannot be ruled out completely.

Should you learn about unwanted side-effects – even if it is doubtful that the side-effect has been caused by our product – please kindly contact us.




Contraindications / interactions

If a patient has known allergies against or hyper-sensitivities towards a component of this product, we recommend not to use it or to do so only under strict medical supervision. In such cases, we will supply the composition of our medical device upon request. The dentist should consider known interactions and crossreactions of the product with other materials already in the patient's mouth before using the product.

**Note**

Please ensure that all above information is readily available to the attending dentist in case this product is used for a special design. Please also note the material safety data sheets available.

Dental ceramics in accordance with DIN EN ISO 6872, type II, class 5.

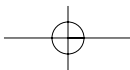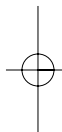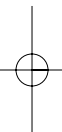CE 0297

SCHÜTZ  **DENTAL**
Micerium Group

Schütz Dental GmbH · Dieselstr. 5-6 · 61191 Rosbach/Germany
Telefon +49 (0) 6003 814-0 · Telefax +49 (0) 6003 814-906
info@schuetz-dental.de · www.schuetz-dental.de